

Product Name:	Cell Scrapers
Product Application:	Specially designed Nest® cell scrape can be used to scrap and collect cells easily and effectively. Push the handle with a slight pressure to change the blade angle. Then push down the handle towards the bottom of the container. Rotate the handle slightly to make the blade twist to right direction.
Product Materials:	HIPS (Blade) + PP (handle)
product Specifications:	Temperature range: stored at room temperature Shelf life: 3 Years after date of production (ensure package is in good)
Sterilization:	Yes. Sterilized by E-beam, Sterility Assurance Level: SAL=10 ⁻⁶ . The product has been irradiated and dosimetrically released based on ANSI/AAMI/ISO 11137
Pyrogens:	Non-Pyrogenic
RNase/DNase:	DNase/RNase free
BSE/TSE:	These products are deemed animal free.
Performance Testing:	Each manufacturing lot is sampled and tested in accordance with standard operating procedures. Appearance inspection: qualified Packaging inspection: qualified
Features	<ol style="list-style-type: none"> 1. Blade is made of HIPS and handle is made of PP 2. Free rotating blade, make the blade twist to right direction 3. Easy-tear sterile packing 4. Sterilized by E-beam, SAL=10⁻⁶ 5. Non-pyrogenic and DNase/RNase free


Product Range

Cat. No.	Spec(mm)		Application	/Case
	Total Length	Blade Width		
710001	220	13	blue handle & transparent blade	100
710011	280	20		

Head Office No. 530, Xida Road, Meicun Industrial Park, Xinwu District, Wuxi, Jiangsu, China
Tel: +86+ 510-6800 6788 Email: info@nest-wuxi.com
Online: www.cell-nest.com

Oversea NEST USA (New Jersey/ Phoenix)
NEST scientific 株式会社 (Yokohama, Japan)
NEST Scientific Europe B.V (Netherlands)
Nest Scientific (MENA) FZE (Sharjah, United Arab Emirates)

Technical Drawing

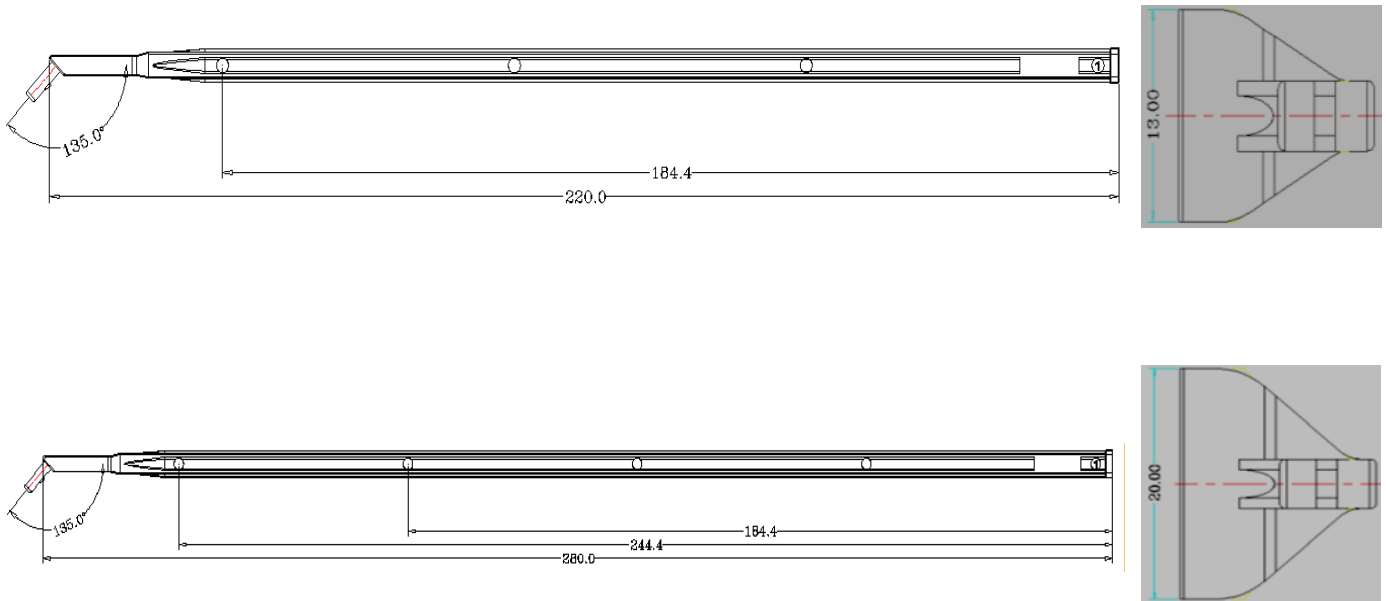

Head Office No. 530, Xida Road, Meicun Industrial Park, Xinwu District, Wuxi, Jiangsu, China
Tel: +86+ 510-6800 6788 Email: info@nest-wuxi.com
Online: www.cell-nest.com

Oversea NEST USA (New Jersey/ Phoenix)
NEST scientific 株式会社 (Yokohama, Japan)
NEST Scientific Europe B.V (Netherlands)
Nest Scientific (MENA) FZE (Sharjah, United Arab Emirates)